

[This package insert must be read carefully prior to use.]

Glucose assay kit  
(Classification No.: 30167002)

## QUALIGENT GLU

### Warning

Measured blood glucose levels may be higher than the actual values in patients with pralidoxime methyl iodide treatment. Therefore, obtain information from the manufacturing authorization holder of pralidoxime methyl iodide concerning its influence on measurement of blood glucose in such patients.

(Measured blood glucose levels may be higher than the actual values in patients with pralidoxime methyl iodide treatment. If hypoglycemic drugs such as insulin are administered on the basis of false high values, serious symptoms of hypoglycemia such as coma may occur.)

### General Precautions

1. This product is for in vitro diagnostic use, and must not be used for any other purposes.
2. Clinicians should make a comprehensive clinical decision based on assay results in conjunction with clinical symptoms and other examination results.
- \* 3. For the effects of an administered drug on the measured value, carefully read the Precautions for Use in the package insert of the drug, especially the section about the effects on laboratory test results. Please also read carefully the “2. Interfering Substances,” in the “Procedural Precautions” section, as well as “2. Precautions for Assessment” in the “Assessment of Assay Results” section, of this package insert.
4. This product should be used only as directed in this package insert. Reliability of results cannot be guaranteed if there are any deviations from the instructions in this package insert.
5. If the reagent accidentally comes in contact with eyes and/or mouth, rinse immediately with ample water as first aid, and consult the doctor if required.
6. Carefully read the operating instructions for each type of automated analyzers prior to using this product. Parameters for each type of analyzers are available, and can be requested from SEKISUI MEDICAL CO., LTD. if required.
7. Perform a quality control test prior to assay to ensure accuracy.

### \* Description (Kit Components)

Component: Ingredients

GLU Enzyme Solution 1:

Glucose-6-phosphate  
dehydrogenase

GLU Enzyme Solution 2:

Hexokinase  
Nicotinamide adenine  
dinucleotide (oxidized form)


### Intended Use

Measurement of glucose in serum, plasma, or urine

### Assay Principle

#### 1. Assay Principle

Glucose is converted to glucose-6-phosphate by the action of hexokinase (HK) in the presence of ATP. Glucose-6-phosphate is converted to 6-phosphogluconic acid by the action of glucose-6-phosphate dehydrogenase (G-6-PDH). At the same time, NAD is converted to NADH, and the absorbance at 340 nm increases. The glucose level is determined by measuring the change in the absorbance of NADH.

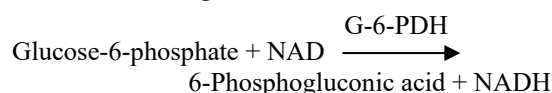

#### 2. Features

- 1) Because this product has the wide measurement range, even urine samples can be measured without dilution.
- 2) Analytical results are hardly affected by interfering substance, because of the rate method.
- 3) This product shows a high sensitivity and reproducibility.

### Procedural Precautions

#### 1. Properties of Samples and Sampling Methods

##### 1) Samples

Serum, plasma and urine may be used.

##### 2) Storage of samples<sup>1)</sup>

The glucose level in whole blood decreases over time. Separate cells from the sample immediately after collection of a blood sample. When serum or plasma samples must be stored after separation, store them at 2-10°C or at -20°C or lower.

Bring samples to room temperature (15-30°C) before use.

Urine samples should be tested on the same day.

#### 2. Interfering substances

Assay results are not affected by free bilirubin (up to 20 mg/dL), conjugated bilirubin (up to 20 mg/dL), hemoglobin (up to 500 mg/dL), formazin turbidity (up to 3000 FTU), or ascorbic acid (up to 50 mg/dL).

#### 3. Others

- \* 1) Always use Serum Multicalibrator (SEKISUI) or Anaserum GLU standard solution for calibration.
- 2) Precautions for assay range  
If the concentration of sample exceeds assay range, dilute the sample with saline and repeat the measurement.

## Dosage/Administration (Assay Procedure)


### 1. Preparation of reagents

Reagent (1): GLU Enzyme Solution 1 is ready to use.

Reagent (2): GLU Enzyme Solution 2 is ready to use.

### \* 2. Assay Procedure

This product is compatible with Hitachi 9000 series and LABOSPECT series automated analyzers. Assay procedure is indicated below.

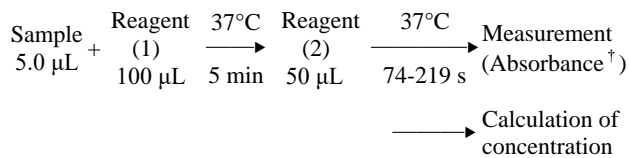

<sup>†</sup>Absorbance: The difference in absorbance between 450 nm and 340 nm

Calibration material: Serum Malticalibrator (SEKISUI) or Anaserum GLU standard solution (Manufacturer's assigned value)

Reagent blank: Purified water or saline

## Assessment of Assay Results

### \* 1. Reference Interval<sup>2)</sup>

Fasting blood glucose level:  
73-109 mg/dL (Within the JCCLS common standard)

### \* 2. Precautions for Assessment

There may be reactions or interfering reactions with non-target substances. If assay results appear to be unreliable, repeat the measurement (if necessary, after dilution) or try another analytical methods.

## Performance

### 1. Sensitivity

- 1) Reagent blank: The change of absorbance is  $\leq 0.005/\text{min}$ .
- 2) Sensitivity: The change of absorbance is 0.030-0.060/min per 300 mg/dL of glucose.

### 2. Accuracy: 90-110% of the expected assay value

### 3. Within-run Reproducibility:

Coefficient of variation  $\leq 3\%$

(Test methods used for 1.-3. are in-house methods.)

### \*\* 4. Measurement Range<sup>7)</sup>: (On a LABSECT 006 automated analyzer)

Serum / Plasma: 1-2000 mg/dL

Urine: 3-7200 mg/dL

### 5. Correlation<sup>7)</sup>

- 1) Serum  $N = 60$   $r = 0.999$   $y = 1.01x - 1.66$   
Control method: Approved in vitro diagnostic (hexokinase method)
- 2) Plasma  $N = 80$   $r = 0.999$   $y = 1.00x + 1.00$   
Control method: Approved in vitro diagnostic (hexokinase method)
- 3) Urine  $N = 50$   $r = 0.999$   $y = 0.99x + 0.29$   
Control method: Approved in vitro diagnostic (hexokinase method)

### 6. Standard Material

SRM917 (NIST)

## Precautions for Use or Handling

### 1. Precautions for Handling (to Ensure Safety)

- 1) All samples used in the test should be handled

as a material possibly infected with HIV, HBV, HCV, or other viruses. To prevent infection, use disposable gloves and avoid mouth pipetting during the test.

- 2) Sodium azide is added as an antiseptic agent in the GLU Enzyme Solution 1 and GLU Enzyme Solution 2. Therefore, if the reagent comes in accidentally contact with eyes, mouth or skin, rinse immediately with ample water as first aid, and consult the doctor if required.

### 2. Precautions for use

- 1) This product should be stored as directed, without freezing. Freezing can deteriorate the reagents, which can produce inaccurate results. Therefore, avoid using the reagents which have been previously frozen.
- 2) Do not use expired reagents. Use of such reagents cannot guarantee the reliability of measurement values.
- \*\* 3) The product should be tightly stoppered and stored as instructed when not in use.
- 4) Do not perform the assay under direct sunlight.

### 3. Precautions for Disposal

- 1) Before disposal, used samples and their containers must be immersed in sodium hypochlorite solution at a concentration of greater than 0.1% for longer than 1 hour or autoclaved at 121°C for 20 minutes.
- 2) To prevent infections from spilled samples or solutions containing samples, wipe the spilled area thoroughly with disinfectants such as sodium hypochlorite solution at a concentration of greater than 0.1%.
- 3) The reagents and treated samples should be discarded as medical waste or industrial waste according to the waste disposal regulations.
- 4) The reagents should be disposed of in accordance with the Water Pollution Control act or related regulations.
- 5) Sodium azide has been added as an antiseptic agent in the GLU Enzyme Solution 1 and GLU Enzyme Solution 2. It can react with lead or copper pipes to produce the highly explosive metal azide. Therefore, the reagent should be flushed with large amounts of water during disposal.

### 4. Other precautions

- 1) Do not use the containers for other purposes.
- 2) Do not take apart the reagent cartridge before.

## Storage and Shelf Life

1. Storage temperature: 2-10°C

2. Shelf life: 10 months from the date of manufacture (The expiration date is printed on the outer package.)

## Packaging

Name	Package		
QUALIGENT GLU	Set (Cassette for Hitachi 9000 series)	GLU Enzyme Solution 1 1 × 15.0 mL	× 2
		GLU Enzyme Solution 2 1 × 7.5 mL	
	L set (Set for Hitachi LABOSPECT series)	GLU Enzyme Solution 1 1 × 40 mL	× 2
		GLU Enzyme Solution 2 1 × 20 mL	

## References

- 1) Sasaki M. et al.: Sampling of constituents of the human body, Tokyo: Kodansha, 1972. Japanese.
- \* 2) Kanai M. (supervising editor): Kanai's manual of clinical laboratory medicine. 35th ed. Tokyo: Kanehara Shuppan. 2020. Japanese.
- 3) Kanai I and Kanai M. (author and editor): Kanai's manual of clinical laboratory medicine. 32nd ed. Tokyo: Kanehara Shuppan. 2005. Japanese.
- 4) Mieling G. E., et al.: Clin. Chem, 25, 1581, 1979.
- 5) Koda K.: Med Tech, 10, 143, 1982. Japanese.
- 6) Nishiyama H. et al.: J Med Pharm Sci, 42(4), 647, 1999. Japanese.
- 7) SEKISUI MEDICAL CO., LTD. In house data. Japanese.

## Contact

SEKISUI MEDICAL CO., LTD.  
international@sekisui.com

## Manufacturer

**SEKISUI MEDICAL CO., LTD.**

**1-3, Nihonbashi 2-chome, Chuo-ku, Tokyo 103-0027 Japan**

"QUALIGENT" and "Anaserum" are trademarks or registered trademarks of SEKISUI MEDICAL CO., LTD. in Japan and/or other countries.

All other brands, product names and service names are trademarks or registered trademarks of their respective companies.